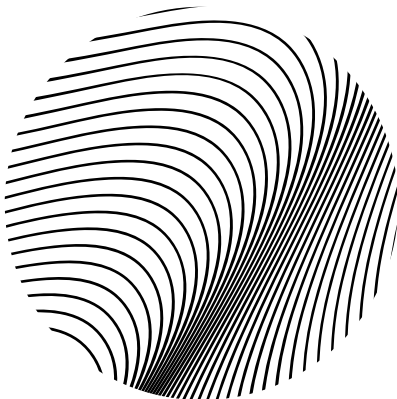
SUPPORTIVE
SHORT
SHORT DE
SOUTIEN

HOW TO WEAR THEM

Prior to using the device, read the following instructions completely.

We're stoked you're stoked to get into your new STOKO SUPPORTIVE SHORTS. Trust us, they're like nothing you've worn before! Here's a couple tips and tricks to make your first experience a little easier.

THEY ARE HARD TO GET ON THE FIRST TIME

Sit down and inch them up like a wetsuit - they will get stuck if you try to pull them on like normal shorts. It could take 5 minutes or more to put them on for the first time, so follow the instructions on the next page which walk you through the essential steps to ensure a flawless fit.

THEY NEED TO BREAK IN

They might feel small at first! Stretch the waistband outward and pull on the cables to distribute support down your thighs. As you move they will mold perfectly to your body, so the more you move, the more the cables will work their magic.

SUPPORTIVE APPAREL FEELS DIFFERENT

Unlike a traditional brace or support, your Supportive Apparel will not feel rigid or bulky. The Stoko Support System is designed to mimic your muscles and ligaments, reinforcing your body and engaging when you need it.

PUTTING ON





1. Pull the dials out to disengage the support system.
2. Stretch the waistband and legs out to relax the cables.
3. Sitting down, pull the shorts over each foot and up past your knees. Position the bottom of the shorts just above the knees.
4. Stand up and pull the waistband over your hips. If the waistband or leg cables feel too tight, pull out on the cables to give yourself some more space. This will be difficult before you have adjusted the support for the first time.
5. Check your alignment: **(I)** The waistband is centered on your body. **(II)** The gusset is centered on the inside of your leg. **(III)** The white Stoko Icon is centered on the outside of your leg.






ADJUSTING

6. Push the dials in and twist to engage support - you will feel them tighten at your waist first.
7. To distribute the support, pull the waistband and thighs outward, then move (squat, jog, etc.).
8. STOKO SUPPORTIVE SHORTS are adjustable and provide more support the more you tighten. Your support is personal; adjust according to comfort or medical recommendation.
9. You will know your shorts are working as you feel it tighten across your waist and thighs.
10. If it feels like the shorts are causing your feet to point outward or inward, reference step 5 to re-align to a neutral position.

REMOVAL

11. Pull the dials out to disengage the support system.
12. Pull the waistband outwards to release the cables and work the shorts down your legs and off.

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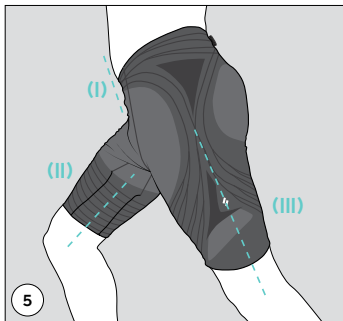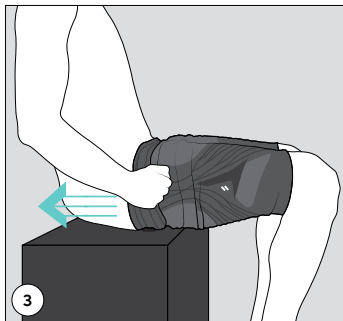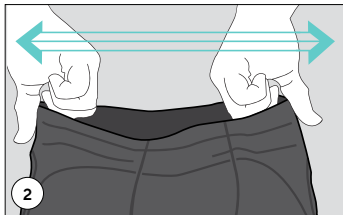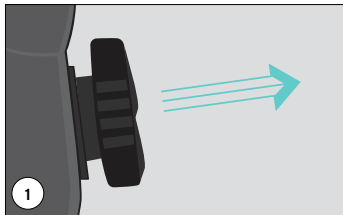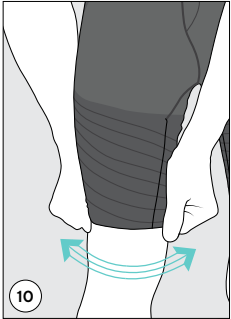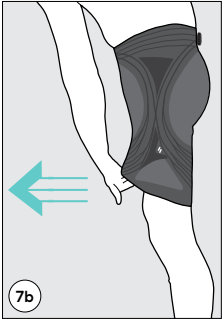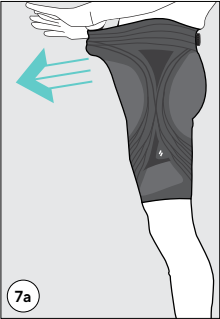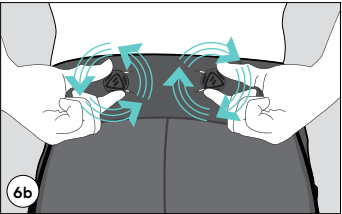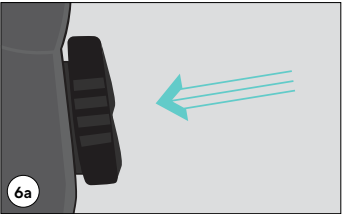THE REST YOU NEED TO KNOW

INTENDED USE

The SUPPORTIVE SHORT is intended for stabilization and support of the lower lumbar, hips, and sacroiliac joint.

INDICATIONS FOR USE

- Lower lumbar pain and strains
- Hip pain, instability, osteoarthritis, and strains
- Sacroiliac (SI) pain, instability, and dysfunction

CONTRAINDICATIONS

If you have any of the following conditions, consult your physician prior to use:

- Pregnancy
- Unstable or displaced fracture of the spine or lower body
- Conditions exacerbated by compression or pressure
- Venous thromboembolism, venous

insufficiency or vascular compromise, impaired blood flow, impaired lymphatic drainage, or soft tissue swelling in the lower body

- Impaired sensation in the spine or lower body

WARNINGS

- Stoko Design Inc. accepts no responsibility for any injury which might occur while using Stoko Design Inc. products or from improper device usage.
- This device will not prevent, reduce or benefit all injuries. Consult your health care provider to ensure the SUPPORTIVE SHORT is suitable for you or if you want to use it together with other products.
- Follow all instructions for use of the SUPPORTIVE SHORT. Failure to properly size, position, tighten and adjust the support will compromise

performance and comfort.

- Do not over-tighten. Any aids applied externally to the body may, if tightened excessively, lead to local pressure points or constrict the underlying blood vessels or nerves.
- In rare cases the STOKO SUPPORTIVE SHORT may cause irritation, discomfort, abrasions, or joint pain. In these cases, discontinue use immediately. If you experience an unusual reaction consult with a healthcare professional.
- Do not use the STOKO SUPPORTIVE SHORT over an open wound.
- Do not modify the product. Stoko Design Inc. accepts no responsibility for the consequences of any adjustments, modifications and/or repairs made to the STOKO SUPPORTIVE SHORT by any user or other third party.
- This product does not support above the lower lumbar spine.
- Remove the SUPPORTIVE SHORT during longer periods of rest. Do not use the support system for longer than 2 hours without a 15-minute rest.

DUTY TO REPORT

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

INTENDED FOR SINGLE PATIENT USE.

USE LIFE

Discontinue device use after 5 years or if you see significant signs of wear or failure.

PRODUCT CARE

Keep your STOKO SUPPORTIVE SHORT fresh! Follow the instructions on the inside of the waistband, wash in the washbag provided, and hang to dry. Do not expose the product to direct heat.

MATERIAL CONTENT

Garment: 64% Nylon / 36% Elastane
Support System: UHMWPE, Nylon, Polycarbonate, Polyoxymethylene, Polyurethane.

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